

# Beautiful Cell Smear



**SPINNER LT1000**

Medical device No: 17B1X10002LP0001



**JAPAN**



# Desktop spreading Spinner LT1000

## ● KEY BENEFITS

- Quick and Easy**      No individual difference and Fixed quality  
 Anyone can prepare the sample of the same quality
- Safe**                      Continuous negative pressure prevents infections
- Good operability**      Desktop type, easy layout.  
 Easy operation equipped with only one button.


## ● HOW TO USE



① the slide glass set




② dropped in the middle of the slide glass (0.25ml blood)

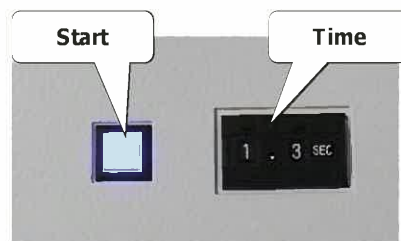

③ set the Time, and push Start button. you can get the cell smear

## ● SPECIFICATIONS

Product name	LT1000
Size	200(W)×251(H)×403(D)mm ※Top cover closed
Weight	12kg
Power supply	100~300VA 50-60Hz
Hardware requirement	15~32℃、30~85%RH
Spinning speed	4700rpm
Spinning time	0.1~9.99 (0.1s)
Option	Foot switch , body color change

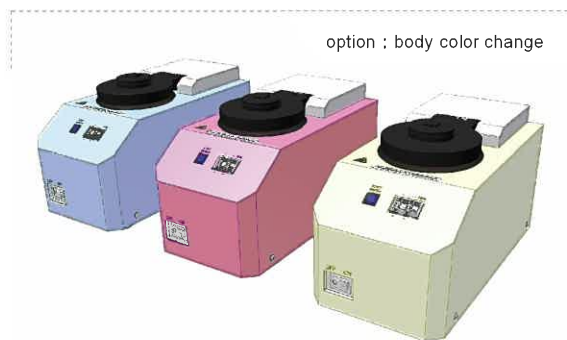

**⚠** When using this device, please read the instruction manual carefully.

### INFORMATION

## LIONPOWER Co.,Ltd

Tsu 5, Tsukidu, Komatsu, Ishikawa,  
923-0972, JAPAN  
TEL +81-761-44-5411 FAX +81-761-44-8080  
E-mail hi@lionpower.co.jp



<http://www.lionpower.co.jp/cart.html>